

# **YF-FR069 Graphene Heating Therapy Waist Belt with Vibration Massage**



## Safety precautions

1. Each using period is 20 minutes, check your waist condition, and then adjust the using times appropriately.
2. It is recommended to use the "low gear" mode at the beginning.
3. Stop using this product right away if you feel pain in your waist or around your waist.
4. If any problems arise in the skin or legs during use, immediately stop using this product and consult a physician.
5. It is recommended not to overuse massage functions.
6. When using this product, your waist may turn red for a short time. This does not mean there are questions. If this symptom persists without relief, consult a surgeon.



Warning Misoperation may lead to a hazard



Caution Wrong operation may lead to rolling or property damage

## Product considerations

This product contains lithium batteries and should not be disassembled, refitted, immersed or thrown in a fire, or it may damage the product.

- Stop using this product immediately if you feel unwell due to its abnormality. Talk to your doctor and follow your doctor's advice if you have unusual sensations in your body or skin.
- Do not store in the reach of children and infants. In addition, do not allow them to touch or handle this product. Otherwise, it may cause burns, electric shock or personal injury.
- Do not damage or modify, or forcefully bend or pull the power cord. Also, do not place heavy objects on or squeeze the power cord.
- Do not use the product in the bathroom or shower.
- Stop using this product right away if you feel pain or discomfort in your skin or other areas of your body.
- Always unplug the adapter from the household power outlet when cleaning.
- Properly dispose of the packing bag of this product away from infants. Do not wear the machine body for any other activities.
- Do not lie or sit on the product.
- It is prohibited to bring great impact and load to products.
- Do not cover quilts or blankets.
- Do not use it while sleeping.

## **Person used with caution**

- Patients with waist nerve disorder or skin sensory disorder.
- Patients injected with metal, silicon or plastic materials into the hand due to fracture or plastic surgery.
- Heat-insensitive people.
- Patients with waist sequela due to traffic accidents, etc.
- Person who is vulnerable to an implanted medical electrical instrument in vivo that is compromised by an electromagnetic malfunction such as a cardiac pacemaker.
- People identified as losing their sense of temperature.
- People with malignant tumors.
- People with heart disease.
- People who are restless regularly at the beginning of pregnancy or shortly after giving birth to a child.
- People with perceptual disturbances caused by diabetes.
- People with skin scars.
- People require quietness and people in poor physical condition.
- People with osteoporosis, severe fractures, sprains, muscle pain, acute (head) pain.
- People under 38°C (fever).
- People under acute inflammatory symptoms, burnout, aversion to cold, blood pressure changes and other symptoms.
- Other people are treated in medical facilities.

## **Using method**

1. Before use, please check whether the product can be started. If it cannot be started, please charge the product with the attached charging line before use.
2. Place the massager over your knee.
3. Press and hold the switch key for 2 seconds, the buzzer drops once, and the default temperature of startup is 40°C.
4. When powering on the product, press the "+" warming key once, add it to 60°C and then press the "+" key, it "beeps" twice (digitally unchanged); press "-" temperature reducing key once to reduce it to 40°C, it "beeps" twice (digitally unchanged), and the temperature range is (40°C - 45°C - 55°C - 60°C).

5. When powering on the product, press M massage mode key, it “beeps” once, the product switches in turn between three massage modes (● blue light - gentle massage, ● green light - soothing massage, ○ white light - strong massage) and off mode, respectively, and the indicator lights correspond to blue - green - white - turnoff - off, respectively, and a warning is dropped for each switching mode buzzer.
6. After powering on, press M key for 2 s to enter screen lock mode, M indicator light is ● red, and it is impossible to switch mode and add or subtract temperature in screen lock state. Press and hold for 2 seconds to unlock (no light display), and a single “beep” will sound each time the screen lock/unlock mode is toggled.
7. In the start-up state, press and hold the power button, the buzzer “beeps” once, and the product will shut down and turn off all the indicators.

## Charging function

When the product is no power, M lamp is ● and flashes three times and “beep” two or three times. When the product is charged, the respiratory lamp ● shows green light, and the ● green light is always on when it is fully charged. If the product needs automatic shutdown during charging, it can not be used while charging.

## Working time

The default time for startup of the product is 20 minutes. When switching off the heating temperature or vibration mode within 20 minutes, it turns off timely and all functions stop, the light off and it shuts down.


## Cleaning and maintenance

### Host care


It is necessary to properly protect the display screen of the host with touch function. For cleaning the host, please use soft cloth to wipe it.

## Product cloth area


Water cleaning is not supported in the cloth area. When cleaning is required, sterilize with alcohol to remove sweat stains or wipe with a clean wet towel.

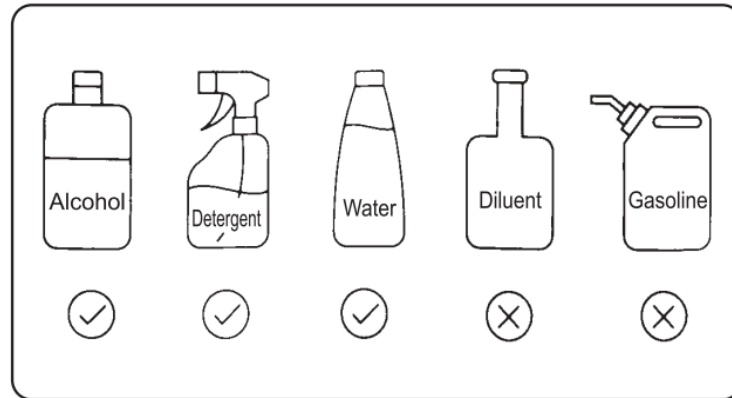

## Key description

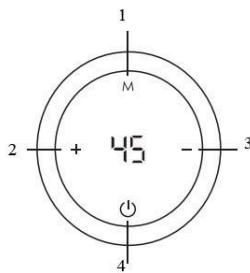

1. Mode key/screen lock
2. Gear
3. Gear
4. Power switch
5. Magnetic suction buckle
6. Magnetic suction buckle
7. Type-C charging port

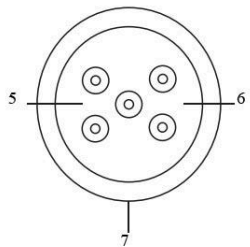

## Failure handling instructions

If the product can't function properly, check that the following condition has occurred. If our advice can't solve the problem, return the product to your distributor for repair.

### **1. Knee protection products cannot be started**

Cause: Product remains connected to a domestic power outlet; battery is depleted; power button is not pressed securely.

Solution: Remove the charging wire from the body; charge; press the power switch key firmly until you press it into place.

### **2. Charging indicator does not on**

Cause: The USB charging cable is not properly plugged into the power outlet.

Solution: Plug the USB charging cable into the power outlet correctly.

### **3. Can only be used for a few minutes even if it has just been charged**

Cause: End of battery life (approx. 5 years)

Solution: Consult your distributor (do not replace the battery on your own)

### **4. Body does not heat**

Cause: This product may not function properly when used outside the recommended ambient temperature range.

Solution: Use within the recommended operating temperature range of 0 ~ 40°C.

### **5. Product stops working during run**

Cause: Battery has run out; after about 15 minutes, the product stops working automatically.

Solution: charge under low battery.

## **Storage requirements**

Keep out of reach of children and away from moisture, heat, direct sunlight and water splashes. Power should be turned off and recharged every two months when not in use for extended periods.

## Use environment

	Use environment	Storage and shipping
Temperature range	0-40C	-10-50C
Humidity range	<80%RH	<80%RH
Atmospheric pressure	70-106KPA	70-110KPA

## Scrap requirements

Batteries after product expiry should be disposed of in accordance with local regulations or environmental regulations.

## Warranty regulations

1. The warranty period for this product may vary from country to country and will be extended as far as permitted by the law of that country.
2. This product is warranted for a period of 1 year from the date of purchase.
3. If quality problems arise during the warranty period, the Company will provide free maintenance or replacement. However, the machine damage caused by your improper use or irregular operation is within the scope of warranty. We will replace your damaged parts with old ones at the ancestral price. Old instruments are owned by the Company.
4. Please show me your product warranty when you need a warranty.

## **Warranty or return may not be offered**

1. Product damage caused by unpredictable events such as fire, earthquake and flood.
2. Instrument malfunction and damage caused by your fault, misuse, use under other abnormal conditions, etc.
3. Instrument damage caused by your self-disassembly and improper storage (such as rat damage, liquid penetration, etc.)
4. Damage caused by forced use of the product in a non-product use environment, including but not limited to exceeding the workload.
5. Damage caused by other circumstances beyond our control.
6. The Company shall not be liable for any incidental damages in connection with the function or use of this Product (such as loss of your business interests, business interruption, other losses caused by the use of instruments, costs caused by the use of this Product requiring replacement of equipment, devices or services, your customer claims, etc.).

The Company shall bear the category and degree of accidents in accordance with relevant laws and regulations.