


CR-250 Smart Pill Organizer with Electronic Timer & 7 Compartments



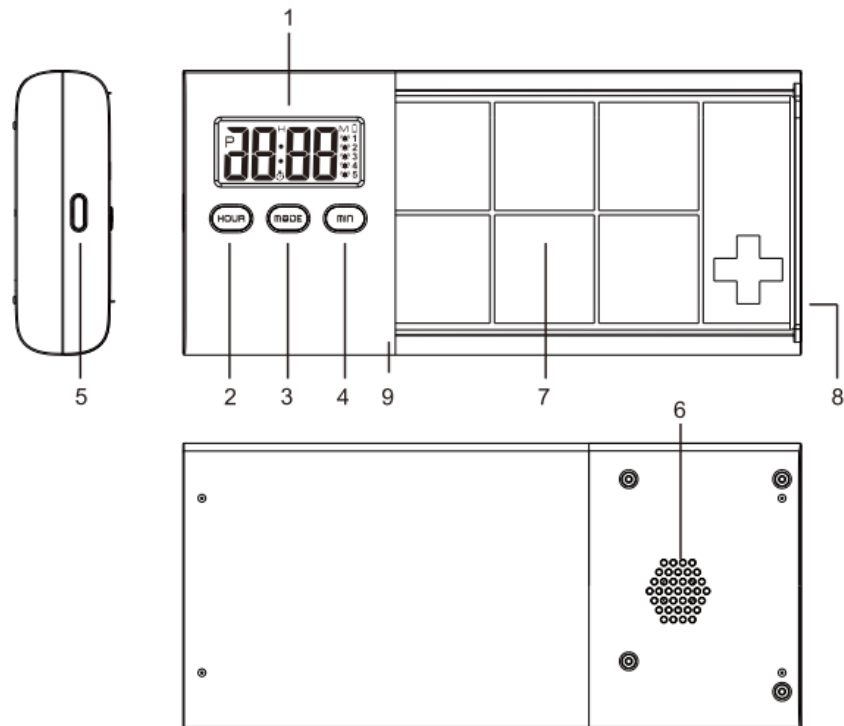

1. Display
2. Hour SET Button
3. Mode / Mode SET Button
4. Minute Set & Start / Stop Button for Timer
5. USB Type-C Charging Port
6. Buzzer Hole
7. Pills compartment
8. Pill Box Ejector Button
9. LED Indicator Light

Features


1. Five Pill Reminder Alarms
2. Alarm Buzzer with Vibration & Flashing Alert Light
3. Seven Pill compartments
4. Simple alarm setting
5. Clock with 12/24 Hours selectable display format.
6. Countdown timer; max. 23 Hours & 59 Minutes.
7. Built-In Rechargeable Battery

Getting Started - First Time Use

Press and hold the MODE button for 2 seconds until the LCD display is activated.

Clock / Alarm / Timer Mode

Press MODE button simultaneously to view the various modes:

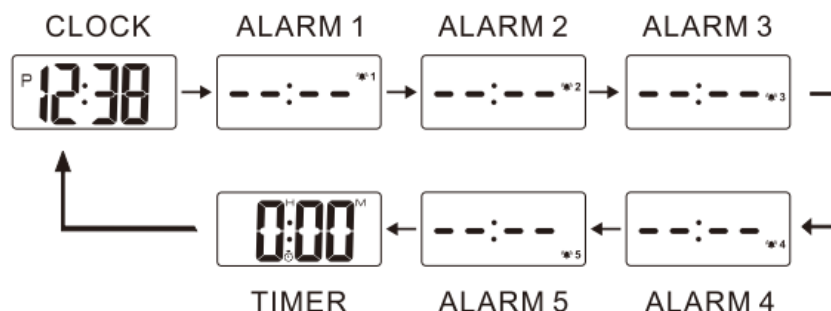

Clock Mode

- Under clock mode, press and hold the MODE button until the digits start blinking 12 Hr / 24 Hr.
- Press Hour or Minute button to toggle select 12Hr or 24 Hr clock format.
- Press the Hour / Minute buttons as many times as necessary to set the current clock time.
Alternatively, press and hold the Hour / Minute button to advance to the set value quicker. Release button when the required value is reached.
- Press the MODE button to acknowledge the set time or leave the unit idle for 1 minute to confirm.

12 Hour clock display



The letter 'P' will show up when the clock time is between 12:00pm to 11:00pm.

Alarm Mode

1. Alarm 1

- Press MODE button until the display shows Alarm 1 icon “🔔1”.
- Press and hold the MODE button until the digits start blinking to enter ‘Alarm 1’ setting mode.
- Press the HOUR / MINUTE buttons to set the alarm time for ‘alarm 1’.
- Once the alarm time is set, press the MODE button to confirm.
- The alarm time 1 will stop blinking. The display will show the alarm time set for alarm 1 and also the icon ‘🔔1’ is visible.

2. To set the alarm for **Alarm 2** (🔔2), Alarm 3 (🔔3), Alarm 4 (🔔4) and Alarm 5 (🔔5), follow the same setting procedures as alarm 1 above.

3. Disable alarm

- Press MODE to select the alarm time from 1 - 5 which you want disabled.
- Press the MINUTE button, the alarm time on the display will change to 4 dash lines (see diagram below). The alarm is now disabled:-



4. Alarm sound / Vibration

When the alarm sound goes off, it will be in the following alarm sound / vibration sequence:

- First 10 seconds: Beep and Vibrate
- Second 30 seconds: Beep at a rapid speed
- Last 20 seconds: Beep and Vibrate at a rapid speed

Press any button to stop the alarm sound and vibration.

Countdown Timer (🕒)

- Press the MODE button simultaneously until the display shows the timer mode.



- Press and hold the MODE button until the digits from the timer start blinking.
- Press HOUR / MINUTE buttons to set the desired countdown time.
- Press the MODE button to confirm the countdown time set.
- Press the MINUTE button to start/stop the countdown timer.
- Once the timer start to countdown, the icon '🕒' will start blinking.
- When the timer reaches its final minute, the display will change and show the countdown timer for the last 59 seconds.



- The timer will beep and vibrate for 1 minute similar to the alarm function.
- Press any button to stop the sound.
- Once the alarm sound is stopped, the display will revert to the previous countdown time set.
- To set a new countdown time, press and hold the MODE button to enter setting mode.

Note: A new countdown time can only be set when the timer has stopped running. If the timer is currently running, press the MINUTE button to stop it.

If you see the timer icon in clock mode, it means the timer is activated and running.




Battery Charging

When the rechargeable battery power is low, the low battery indicator at the top right hand side of the display will show up. It is time to recharge the battery.



- Plug one end of the charging cable to the unit's type-C socket and the other end to an USB power adaptor or computer USB port.
- The blue indicator light will turn ON to inform the user the unit is charging:

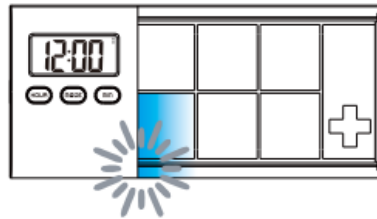

- Once the unit is fully charged, the blue light will turn OFF.

Pill Box

Press the Pill Box Ejector button (the side panel of the pill box tray) to eject the pill box tray out from the unit.



Maintenance

If the unit is not in use for more than 6 months, it is recommended to charge the unit every 6 months to prolong the battery life.